LANSINOH®

Affinity™

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Alexandria, VA 23314

1.800.292.4794

Lansinoh® is a registered trademark of Lansinoh Laboratories, Inc.

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DE52015IF1009

Quick Start Instructions

If you are using the Lansinoh Affinity Breast Pump for the first time, it is important to read the full instruction manual before you begin.

Quick-Start Instructions

Start with Phase I (Let-Down Phase)

Turn on power by pressing the on/off button 

This will take you into Phase I – the Let-Down Phase

Phase II (Expression Phase)

The pump will automatically switch into Phase II after two minutes.

If you want to bypass or reduce Phase I, press the Let-Down button. 

The pump will switch to Phase II immediately.

Adjust suction

While you are in Phase I or Phase II, you may adjust suction level by pressing (+) or (-) next to the suction icon. 

Adjust speed

While you are in Phase II, you may adjust speed level by pressing (+) or (-) next to the suction icon. 

Change to Phase I (Let-Down Phase)

At any time during the Expression Phase, you can switch back to Phase I (Let-Down Phase) by pressing the Let-Down Button. 

Stop pumping

Turn the pump off with the on/off button. 

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Nivel 1 (L1) a Nivel 6 (L6)

www.lansinoh.com

Congratulations On Your New Lansinoh Affinity Double Electric Breast Pump!

The Affinity Breast Pump provides many of the features moms are looking for in a high-quality double electric Breast Pump.

- Allows for single or double pumping to save time
- Lansinoh Custom Expression™ technology allows mom to separately adjust the speed and suction settings for comfortable and efficient milk expression
- Provides two phases of expression – the Let-Down Phase and the Expression Phase – to simulate the natural rhythm of breastfeeding
- Guaranteed to prevent milk back up in tubing and pump
- AC or battery operated—two options to fit your needs
- Small, quiet and lightweight

If you have any questions, please call our customer service department at **1.800.292.4794** or email **customerservice@lansinoh.com**.

Contents

Quick Start Instructions	1	Pumping Tips	10
Important Safeguards	4	Using The Affinity Breast Pump	10
Pump Diagram	5	Power Options	12
Affinity Breast Pump Parts	7	Breastmilk Storage Guidelines	13
Cleaning Instructions	8	Troubleshooting	14
Assembling The Breastmilk Expression Set	9	Technical Data	14
		Warranty	15

IMPORTANT SAFEGUARDS

The Affinity Breast Pump is a personal care item and for hygienic reasons should never be shared between mothers or resold from one mother to another. Either of these will void the warranty.

If you are a mother who is infected with Hepatitis B, Hepatitis C or Human Immunodeficiency Virus (HIV), pumping breastmilk using the Lansinoh Affinity Breast Pump will not reduce or remove the risk of transmitting the virus to your baby through your breastmilk.

This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to other devices in the vicinity. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to other devices, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- ▶ Reorient or relocate the receiving device.
- ▶ Increase the separation between the equipment.
- ▶ Connect the equipment into an outlet on a circuit different from that to which the other device(s) are connected.
- ▶ Consult the manufacturer or field service technician for help.

When using electrical products, especially when children are present, basic safety precautions should always be maintained, including the following:

READ ALL INSTRUCTIONS BEFORE USING

Danger: To reduce the risk of electrocution:

- ▶ Always unplug this product immediately after using.
- ▶ Do not use while bathing.
- ▶ Do not place or store this product where it can fall or be pulled into a tub or sink.
- ▶ Do not place or drop pump into water or other liquid.
- ▶ Do not reach for an electrical product that has fallen into water. Unplug immediately.


Warning: To reduce the risk of burns, electrocution, fire, or injury to persons:

- ▶ This product should never be left unattended when plugged in.
- ▶ Close supervision is necessary when this product is used near children or invalids.
- ▶ Use this product only for its intended use as described in this manual. Do not use attachments not recommended by the manufacturer.
- ▶ Never operate this product if it has a damaged cord or plug, if it is not working properly, if it has been dropped or damaged, or dropped into water.
- ▶ Keep the cord away from heated surfaces.
- ▶ Never use while sleeping or drowsy.
- ▶ Never drop or insert any object into any opening or tubing.
- ▶ Do not use outdoors with the AC adapter or operate where aerosol spray products are being used or where oxygen is being administered.

SAVE THESE INSTRUCTIONS
along with your receipt for future reference.


Pump Diagram

Basic Controls, Display And Features

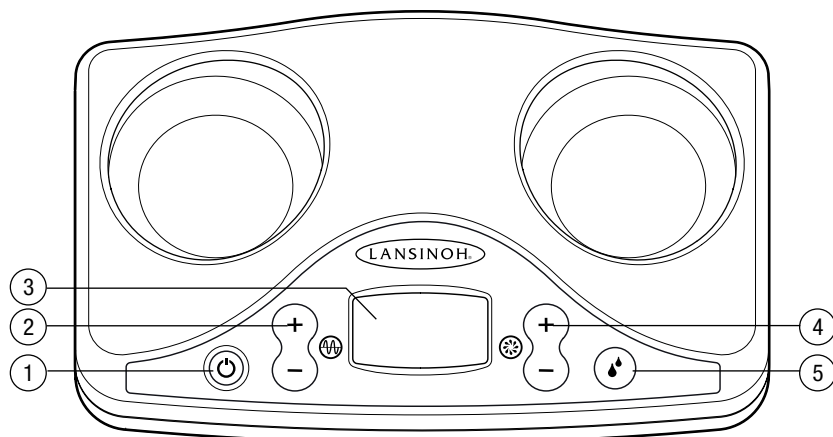
- 1. Power, On/Off button
- 2. Speed increase/decrease button
- 3. LCD screen
- 4. Suction increase/decrease button
- 5. Let-down button

- Power, On/Off Icon
- Cycle/Speed icon
- Suction/Vacuum icon
- Let-down icon

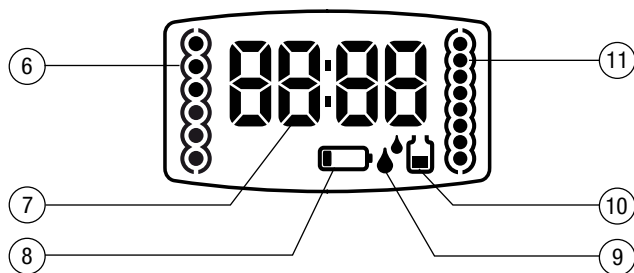
- 6. Speed indicator
- 7. Time/level indicator
- 8. Low battery indicator
- 9. Let-down/Stimulation Phase indicator
- 10. Expression Phase indicator
- 11. Suction indicator

Tubing Adapter Input

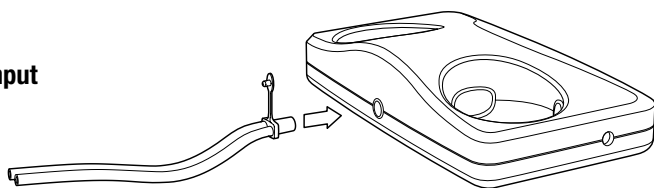
Pump Diagram

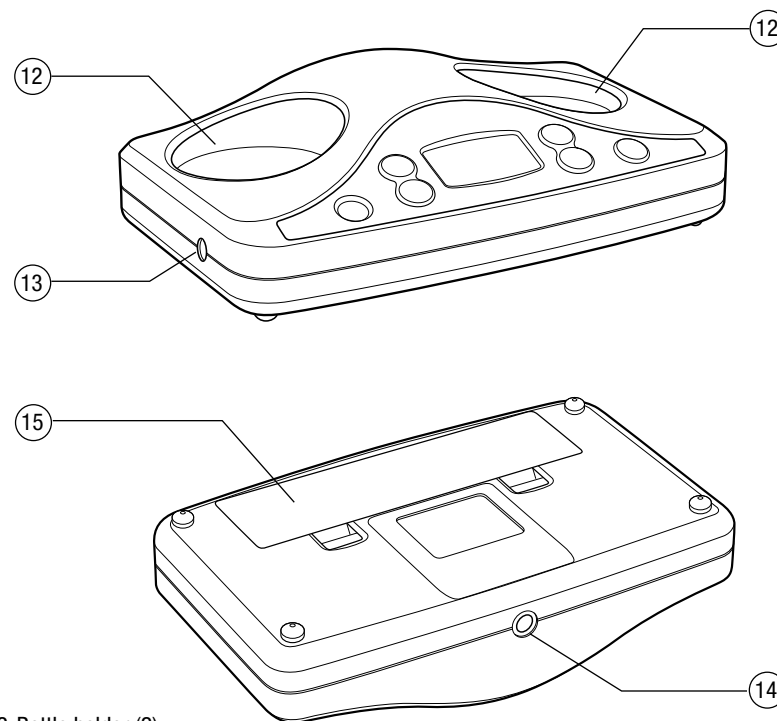
- 12. Bottle holder (2)
- 13. AC adapter input
- 14. Tubing adapter input
- 15. Battery compartment for 6 AA batteries (not included)

AC Adapter

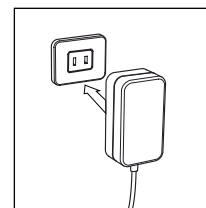
AA Battery Compartment (underneath pump)

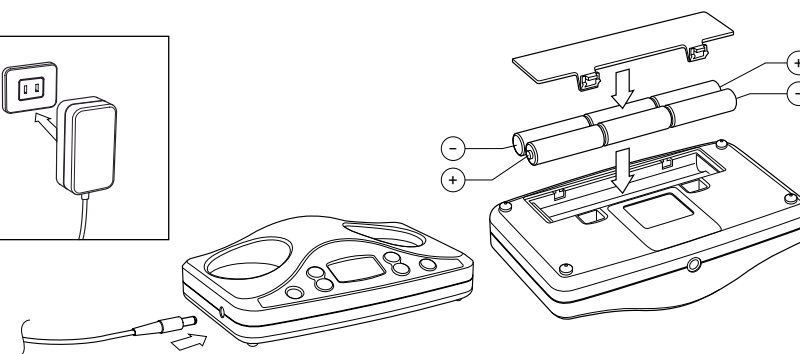Affinity Breast Pump Parts:

Item:	Item Number:
a. Tubing	#50403
b. Breastmilk Expression Set includes:	
b1. ComfortFit™ Breastflange (2)	#50408
b2. Breastflange Body (2)	#50415
b3. Diaphragm (2)	#50402
b4. Diaphragm Cap (2)	#50410
b5. White Valve (2 + 2 extra)	#50401
b6. 5 oz (150ml) Bottles (2)	#50414
c. Breastmilk Storage Bottle Ring (2) and Sealing Disk (2)	#50406 & #50413
d. AC Adapter/Power Cord	#50409
e. Velcro Strap for Tubing	#50417

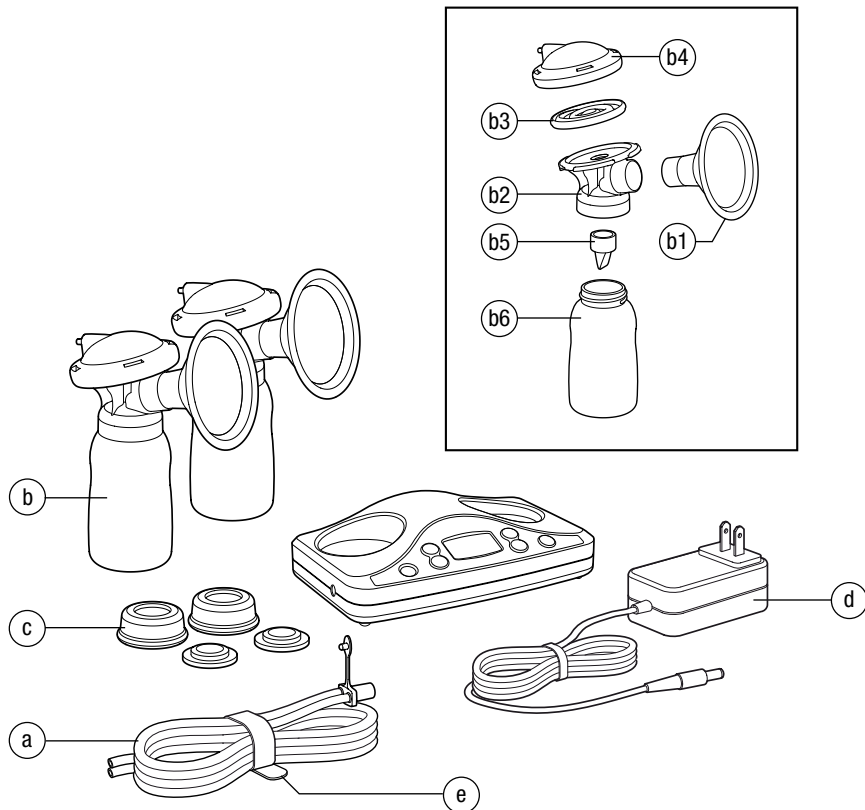Two extra valves are included. Breast Pump carry bag, larger size flanges, cooler and cooler pack are separately available for sale. To order these items or replacement parts, please visit www.lansinohstore.com or call 1.800.292.4794 and reference the item numbers above.

Before Each Use...

- Always wash hands thoroughly with soap and water before handling the Breastmilk Expression Set, the pump, or other attachments.
- Take the Breastmilk Expression Set apart to clean.
- Ensure that the tubing is free of any liquid. Do not operate the pump if liquid is visible in the tubing.

Cleaning Instructions

Unless otherwise instructed by your healthcare provider, follow these sanitizing and cleaning instructions.

Sanitization Prior To First Use

1. Prior to using the set for the first time, disassemble and boil all parts of the Breastmilk Expression Set for 10 minutes.
 - **Note:** Ensure you use enough water to keep parts afloat while boiling.
 - **Note:** Do not boil the tubing and tubing connector.
2. Remove all parts from water immediately after boiling.
3. Completely air or towel dry before assembly and next use.
 - **Note:** Partially dried parts may affect operation of the pump. Please ensure parts are completely dry before use.

Cleaning


- Set aside the tubing and tubing connector. These items do not need to be washed or sanitized because the diaphragm prevents them from coming into contact with the breastmilk.
- All parts of the Breastmilk Expression Set may be cleaned in the upper compartment of your dishwasher or hand washed with dish soap in hot water. DO NOT use solvents or abrasives. Rinse with hot, clean water.
- To clean your Breast Pump, wipe the base with a damp cloth, using only water or mild soap.
- **Do not wash the tubing or tubing connector.**

Note: Before using the Breastmilk Expression Set, visually inspect the individual components for cracks, chips, tears, discoloration or deterioration, as they may affect the functionality of the Breastmilk Expression Set. Please visit www.lansinohstore.com or call 1.800.292.4794 for further guidance and replacement parts.

Assembling The Breastmilk Expression Set

Your Breastmilk Expression Set comes preassembled in Lansinoh product packaging.

Clean Breastmilk Expression Set parts (except tubing and tubing connector) as described in the “Cleaning Instructions” section (see page 6). Once all parts are clean and dry, assemble the expression set as follows:

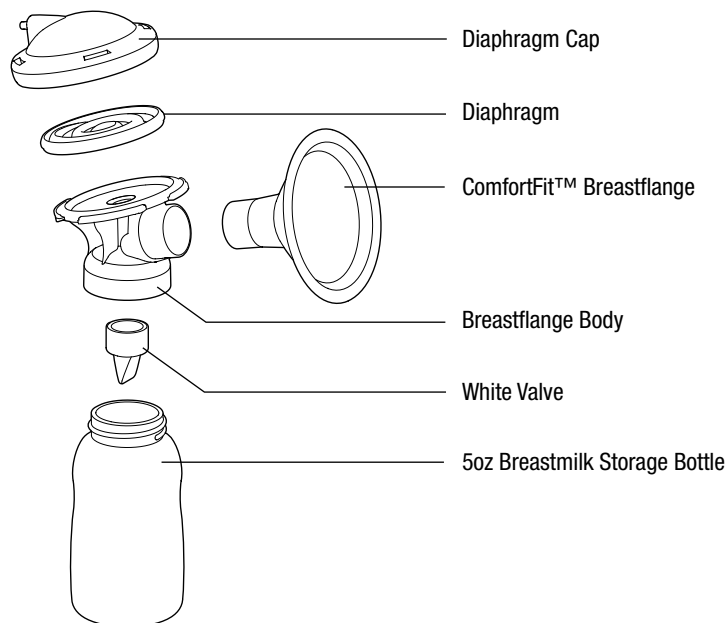1. Place purple diaphragm on top of breastflange body, twist clear diaphragm cap onto breastflange body, enclosing the purple diaphragm. Diaphragm can be used on either side. Ensure diaphragm cap is twisted tightly into place.

2. Push breastflange into breastflange body.


3. Gently secure white valve into the bottom of breastflange.

Note: Use caution when handling and cleaning the small white valve (see page 6).

4. Screw breastmilk storage bottle onto breastflange body.

5. Insert tubing into back of clear diaphragm cap, ensuring a tight fit.

6. Insert tubing connector, at end of tubing, into tubing adapter input on the pump.

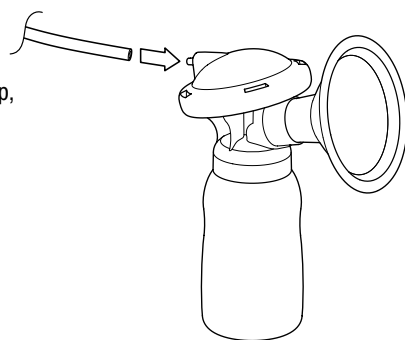Pumping Tips

The most important part of successful pumping is the let-down reflex. This is when your milk begins to flow freely. Warm wet compresses applied to the breasts and light, circular massaging of the breasts prior to and/or during a pumping session will help this natural reflex occur. Thinking of or looking at your baby, or a photo of your baby, may also be helpful.


Some moms can achieve more than one let-down during a pumping session. You may be able to achieve a second or third let-down during your pumping session by pressing the Let-Down button.

How Long To Pump


Double pumping, compared to single pumping, will reduce the time you spend pumping. Remember, successful pumping is a learned skill. Breastmilk is produced on a supply and demand basis; expressing your breastmilk regularly stimulates and maintains your milk supply. If breastfeeding problems occur, call your physician, an International Board Certified Lactation Consultant, or a La Leche League Leader.

Using The Lansinoh Affinity Breast Pump

For Double Pumping

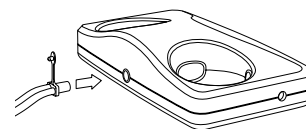
1. Make sure both tubes are attached to the tubing connector and insert into pump base.

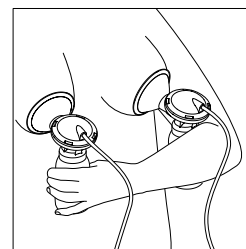2. Hold the Breastmilk Expression Set to your breasts so that your nipples are centered

3. Press the On/Off button to begin pumping. Lansinoh Affinity Pump will start in the Let-Down Phase

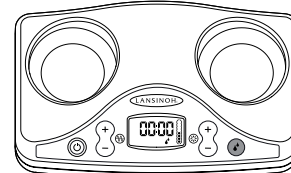If you prefer to bypass the pre-set Let-Down Phase or if your milk lets down sooner than two minutes, press the Let-Down button


Note – the Expression Icon will appear in the display to signify that you are in the Expression Phase.

Using The Affinity Breast Pump

- Adjust the suction to your comfort level while pumping by pressing the increase (+) or decrease (-) button next to the suction icon.

After two minutes, the Let-Down Phase will automatically change to the Expression Phase, indicated by the Expression Icon on the display.

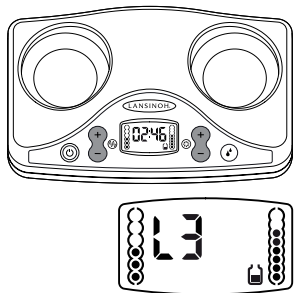
- While in the Expression Phase, you can separately adjust both the speed and suction settings. Choose the setting that is most comfortable for you and results in the greatest milk expression.

Note: The speed or suction level will appear on the display panel ["L1" indicates the lowest level and "L6" (speed) or "L8" (suction) indicates the highest level.] This display can help you remember what speed and suction level is right for you.

- When a pumping session is finished, turn the pump off with the on/off button. After turning the pump off, the display will remain on for a few seconds. Make sure there is no suction on the nipple(s) before removing the breastflanges. If necessary, break the suction by inserting a finger between the breastflange and the breast.

Note: LCD display light will automatically switch off if no buttons are pressed for 20 seconds in order to save power. If the pump is left on but unattended for more than 30 minutes, the pump will automatically power off.

For Single Pumping:

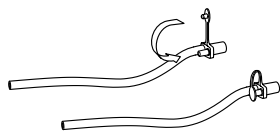Remove one tube and close the tubing connector. Make sure that the plug is placed into the pipe outlet properly and pushed in securely.

Using The Affinity Breast Pump / Power Options

The Adjustable Difference

Your Lansinoh Affinity Breast Pump allows you to adjust both speed and suction settings. Choose the settings that are most comfortable for you and result in the greatest milk expression.

Speed Buttons

These control the speed of the pump. Push the + button to increase the speed. Some mothers have found that their milk flows more quickly if the pump speed is initially set on "high." Push the (-) button to decrease the speed to the highest setting comfortable for you. Find the speed that most closely corresponds to the sucking rhythm of your baby. Your pump can be adjusted from speed levels 1 to 6.

Note: Some mothers find it helpful to begin pumping on a high setting, then lower the setting once milk starts to flow.




Suction Buttons

These control the suction strength. Push the (+) button to increase the suction. Continue to push the (+) button to increase the suction to the highest level comfortable for you. Push the (-) button to decrease the suction. Your pump can be adjusted from suction levels 1 to 8.

Note: If the suction on your breast reaches uncomfortable levels, reduce the suction level by pushing the (-) Suction button, or disconnect the pump connector from its receptacle on the electric pump.

Power Options

The pump may be operated under any of the following power supply modes:

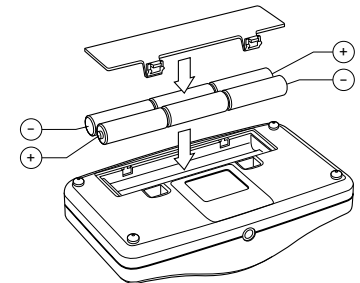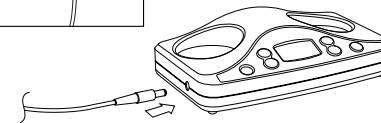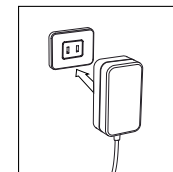AC Adapter (included)

Connect the AC adapter to the power receptacle located on the left side of the pump. Plug the AC adapter into an available electrical outlet. Electrical safety can only be achieved when the Breast Pump is connected to the AC Power Adapter supplied with the product. (See "Technical Data" section, page 12.)

6 AA Alkaline Or Rechargeable Batteries (not included)

Open the battery cover under the pump and insert 6 AA batteries. Make sure the polarity (+) and (-) is correct as indicated. Close the compartment by gently pressing the cover until it snaps into place. Remove batteries if pump is not in use for a long period of time. Dispose of batteries according to regulations for environmental protection and recycling. Replace batteries as a simultaneous set.

Breastmilk Storage Guidelines

Power Options (continued)

If batteries are in the pump, but AC adapter is in use, the pump will use the adapter to power the pump, not the batteries. For the most economical use, we recommend using your pump with the AC adapter provided. Your AC adaptor has a universal AC input and can be used with 120Vac and 230Vac systems with 50Hz and 60Hz mains. The Affinity pump is intended to only be used with the provided external adapter (Lansinoh part # 50409) by BridgePower Corporation, Model JMW110-KA-09-00-B-07, with electrical rating 100-240 V AC, 50-60Hz.

Breastmilk Storage And Thawing Guidelines For Healthy Term Babies

Check with your healthcare provider or an International Board Certified Lactation Consultant for specific storage instructions. When you are freezing your breastmilk, do not fill containers more than ¾ full to allow for expansion from freezing. Date containers with the date you expressed your milk for future reference.

	Deep Freeze	Refrigerator Style Freezer (Self-contained freezer) approx. 0°F	Refrigerator 32°F to 39°F	Room Temperature 66° to 78°F
Fresh Breastmilk	6+ months	6 months	Up to 8 days	Up to 6 hours
Thawed Breastmilk	Do not refreeze	Do not refreeze	24 hours	1 hour

Source: La Leche League International, 2008

Note: DO NOT store breastmilk in freezer door. The door tends to be the warmest part of the freezer.

Date all breastmilk and use the oldest breastmilk first.

If you purchased the optional Cooler Pack and Freezer Pack as accessories to your Affinity Breast Pump, you may safely store up to six (6), 5 oz. bottles in the room at an ambient room temperature of 66° to 78°F for up to ten (10) hours. Freezer packs should be removed once thawed and refrozen.

Thawing

Thaw frozen milk in the refrigerator or place the containers under cool running water, gradually warming the water until the milk is thawed. Thawed milk can be kept in a refrigerator for up to 24 hours.

NOTE: DO NOT refreeze milk once it is thawed.

Gently shake thawed milk before feeding to baby to mix the layers that have separated. Breastmilk is not homogenized and the “cream” will rise to the top of the bottle. Breastmilk can vary in color and consistency depending upon the time of the day that the milk was expressed, what you have eaten, and the age of the baby at the time of pumping.

NOTE: DO NOT MICROWAVE HUMAN MILK! Microwaving alters the composition of the milk. Microwaving also may create “hot spots” in the milk that could burn the infant.

NOTE: DO NOT thaw breastmilk in hot or boiling water.

Troubleshooting / Technical Data

Troubleshooting The Lansinoh Affinity Breast Pump

Interruption of Electrical Service

If you experience an interruption in the normal operation of the electrical and/or mechanical functioning of the pump, simply do the following:

- Unplug the AC adapter cord from the electrical outlet.
- Wait a few seconds, and then plug the AC adapter cord back into the electrical outlet.

Discomfort While Pumping

If you feel any discomfort, you may need a different size flange. Visit www.lansinohstore.com or call 1.800.292.4794.

Decreased Suction

If it feels like the suction has decreased, try the following:

- Remove the tubing and tubing connector from the back of the pump and place a finger over the pump outlet. If you feel suction, the pump is operating properly, but the Breastmilk Expression Set may not be correctly assembled or the white valve may be damaged. If no suction is felt on your finger, try unplugging the adapter from the electrical outlet, wait a few seconds, and plug it back into the electrical outlet.
- Make sure the Breastmilk Expression Set has been properly assembled. Check the assembly information in the section “Assembly Instructions”, page 7.
- Inspect the white valve for a tear or pinhole. This component is critical to achieving proper suction. Replace the white valve with one of the spare valves provided with the product.
- When single pumping, check the tubing adapter on the back of the pump to make sure it is closed for single pumping. If open, close the adapter.
- If double pumping, ensure that the breastmilk expression set is on both breasts before starting the pump.

Increased suction level for single pumping:

If it feels like the suction is higher than usual when single pumping, check that the plug for the pipe is clean and that the airway in the plug is not blocked.

Breastmilk Build-up in Breastflange or Breastflange Body

If breastmilk does not appear to drain fast enough into the bottle, reduce the suction setting on the pump to allow the milk to drain. Once milk has drained into the bottle, you may gradually increase the suction level.

Technical Data

Power Requirements

AC Adapter (included)
AC Adapter – North American Plug
Configuration
110Vac/220Vac wall plug type AC adaptor,
galvanically isolated adaptor with
9Vdc output voltage
Batteries (not included)
6 AA batteries, 1.5V alkaline or
1.2V rechargeable

Regulatory Compliance

UL and CSA approved with North
American power adapter.
UL File Number E323214
CSA File Number E323214
MEDICAL EQUIPMENT WITH RESPECT TO ELECTRICAL
SHOCK, FIRE AND MECHANICAL HAZARDS ONLY IN
ACCORDANCE WITH 60601-1 and CSA C22.2 No.
60950-1 CONTROL NUMBER: 44TS

Type BF Applied Part

Limited Consumer Product Warranty (United States)

This product is warranted to the original consumer purchaser (the "purchaser") to be free from defects in material and workmanship which are not commercially acceptable for the period of one year from the date of purchase on the pump mechanism and the period of 90 days on all other product components. Warranty coverage terminates if you sell or otherwise transfer this product to another person. This warranty gives you specific legal rights and you may also have other rights, which vary, by location.

LANSINOH MAKES NO EXPRESS WARRANTY OF ANY KIND REGARDING THE PRODUCT OTHER THAN THOSE WARRANTIES EXPLICITLY SET FORTH HEREIN. ANY IMPLIED WARRANTY, INCLUDING ANY IMPLIED WARRANTY OF MERCHANTABILITY OR ANY IMPLIED WARRANTY OF FITNESS FOR A PARTICULAR PURPOSE, TO THE EXTENT PERMITTED BY LAW, SHALL BE LIMITED IN DURATION TO A PERIOD OF 90 DAYS FROM THE DATE OF PURCHASE BY THE ORIGINAL PURCHASER.

In the event the product is found by Lansinoh to not meet the above limited warranty, as purchaser's sole and exclusive remedy Lansinoh will repair or, at the option of Lansinoh, replace this product without charge for such replacement or parts or labor. The purchaser shall bear all expense for returning this product to Lansinoh. This warranty does not apply to any product that has been subjected to misuse, abuse, or alteration. Only Lansinoh brand replacement parts are recommended for use with Lansinoh breastfeeding products. Improper or incorrectly performed maintenance or repair, or use of non-Lansinoh brand replacement parts may invalidate this warranty.

TO THE EXTENT PERMITTED BY LAW, LANSINOH SHALL NOT BE LIABLE FOR ANY INCIDENTAL OR CONSEQUENTIAL DAMAGES INCLUDING, BUT NOT LIMITED TO, REPLACEMENT COSTS RESULTING FROM THE BREACH OF ANY WRITTEN OR IMPLIED WARRANTY.

If you wish to make a claim under this warranty, you must return the product to Lansinoh Laboratories Inc., prepaid, together with your dated bill of sale or other dated proof of purchase and a brief statement of the problem, along with a Return Authorization (RA) number. You may obtain the appropriate Return Authorization number and shipping address by calling the phone number or sending an email to the email address provided below.

Items returned without an associated Return Authorization Number will be returned to sender.

To obtain the address to Return goods to: 1.800.292.4794
customerservice@lansinoh.com

To obtain a Return Authorization Number: 1.800.292.4794
customerservice@lansinoh.com

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©2009 Lansinoh Laboratories, Inc. All Rights Reserved. / Todos los derechos reservados.